

iScale[®] Plus

Cellular-Connected Body Weight Measuring System

One Solution for Pediatrics to Bariatrics

Remote Patient Monitoring (RPM)

Manage your most critical chronic care patients using real-time, immediately transmitted weight data.

Patient Friendly, Easy to Use, Automatic Data Upload

No Bluetooth®, no synching, no apps required.
Helps reduce user errors and delays to providing care.

Connected Care Between Visits

Improved patient care, adherence and outcomes.

Improve Clinical Efficiency

Effectively monitor patients remotely using the iScale body weight measuring system.



Connected by



Features

- Weight Capacity Up to 551 lbs
- Larger Footprint
- Large Display Screen
- Non-Slip Surface
- Easy to Clean
- Passive Reading Storage
- Suitable for Ages 10+
- Enhanced cellular connectivity to our Private Data Network via multicarrier SIM technology, connecting to the strongest cellular network available in the device area
- 4 x AA Batteries
- 100% Quality Guarantee

www.smartmeterrpm.com | 1-813-773-4080



SMARTMETER
Engaging Health Data



iScale[®] Plus

Specifications Summary



iScale [®] Plus Product Features and Information		
Display	Digital LCD	
Battery	Four AA (1.5V) batteries	
Dimensions	14 in x 14 in x 1.5 in (360 mm x 360 mm x 37.2 mm)	
Data Transmission	Cellular connectivity to our Private Data Network via multicarrier SIM technology	
Scale Weight	6.7 lb (3.04 kg)	
Measurement Unit	Pound / Kilogram / Stone	
Measurement Range	11 lb to 551.2 lb / 5 kg to 250 kg / 0 st:11 lb tp 39 st: 5 lb	
Accuracy	-Weight Measurement (lb)	-Weight Measurement (kg)
	110 lb: +/- 0.5 lb	50 kg: +/- 0.2 kg
	220 lb: +/- 0.7 lb	100 kg: +/- 0.3 kg
	330 lb: +/- 0.9 lb	150 kg: +/- 0.4 kg
	440 lb: +/- 1.1 lb	200 kg: +/- 0.5 kg
	551 lb: +/- 1.8 lb	250 kb: +/- 0.8 kg
Operating Conditions	Temperature: 41°F to 104°F (5°C to 40°C) Relative humidity: <=90%, Atmospheric pressure: 860 hPa to 1060 hPa	
Storage Conditions	Temperature: -4°F to 140°F (-20°C to 60°C) Relative humidity: <=90% Atmospheric pressure: 500 hPa to 1060 hPa	
Power On	Automatic, "SENSE ON" technology	
Power Off	Automatic, after 10 seconds while showing 0.0	
Mode of Operation	Continuous operation	
Water Protection	IP22 – protected against solid foreign objects >=12.5mm and vertically falling water droplets when enclosure title <=15°	
Degree of Protection	Type BF applied part (top surface)	
Scale Placement	Flat, hard surface	
Measurement Timing	At least 2 hours after waking up or dining	
Use Location	In home	
Age of Use	Adults, Children (10-17 years old)	
Device Lifetime ¹	Five (5) years	
Warranty	Two (2) years	

¹ Device lifetime is not the same as the warranty period. Please consult product agreement for warranty details.



iScale[®] Plus

Specifications Summary



iScale [®] Plus Display Icons	
kg	Weight measurement unit – kilogram
st lb	Weight measurement unit – stone
lb	Weight measurement unit – pound
	Signal symbol – no signal
	Signal symbol – strong
	Signal symbol – medium
	Signal symbol – general
	Signal symbol – weak
	Data transmitting
	Battery level

iScale [®] Plus Display Errors	
	Overload error – retake measurement, if issues persist check if your weight is within range (11 lb to 551.2 lb), stop using this scale for measurement if not in range
	Low battery error – replace all four AA batteries at same time
	Failure to pair-up with network, data transmission error – move scale to an area with a strong signal
	Server API call was unsuccessful – contact customer support
	Query signal timeout – reinstall the batteries, move scale to an area with a strong signal
	No SIM card information was found – contact customer support
	Server API replied certificate is abnormal or device does not exist (certificate error) – contact customer support
	Eeprom checks failed (hardware error) – contact customer support

Customer Support – Call 844.445.8267 or email support@iglucose.com

Why Smart Meter is the No. 1 Supplier of Cellular RPM™ Solutions

The Only Trusted Private Data Network for RPM in the U.S.

- Smart Meter's Private Data Network for RPM is an exclusive ecosystem with robust security layers protecting patient-generated health data from our proprietary monitoring devices to our SmartPartners™ platforms.

Cellular RPM™ is Better Than Bluetooth®

- Patients have problems using Bluetooth – it's that simple.
- Our proprietary cellular-enabled devices deliver the most patient-friendly experience available supporting more frequent testing, additional revenue, an improved volume of health data facilitating care management, and higher customer and patient satisfaction.

Everything Working Seamlessly in the Background

- We simplify the deployment and management of remote patient monitoring so you can focus on transforming healthcare.
- SmartPartners have exclusive access to curated monitoring devices, real-time RPM data stacks, shipping and logistics management, automated refills, and success training services.

100% Quality Guarantee

- We stand behind our solutions: products, data, platform, and services.
- If anything fails to meet your standards, we will resolve the problem to your satisfaction. Guaranteed!


Grow With Us

Whether you're implementing an RPM system for the first time or need to reach new customers with expanded solutions, we're here to help.

www.smartmeterrpm.com | 1-813-773-4080



70% of people have trouble sending health data with BlueTooth® devices.

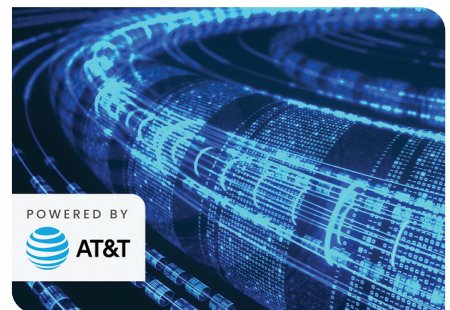

SMARTMETER
Engaging Health Data

5501 W. Waters Ave.
Suite 401
Tampa, FL 33634
www.smartmeterrpm.com
813-773-4080
SC-S-iSCPRODOV-U-0924